



 **Case Medical**  
*Sustainable Products for Instrument Processing*





# Table of Contents

<b>About Case Medical</b>	
Company Overview	2
What We Offer	3
Our Vision	4
<b>SteriTite®   Universal Sealed Containers</b>	<b>5</b>
Full Size Sealed Containers	6
Three-Quarter Size Containers	7
Half Size Containers	8
Long & Wide Containers	9
Mini, Narrow & Narrow Long Containers	10
FlashTite® Containers & Baskets	11
<b>MediTray®   Instrument Baskets</b>	<b>12</b>
Full Size Baskets	13
Three-Quarter Size Baskets	14
Half Size Baskets	15
Long & Wide Baskets	16
Mini, Narrow & Narrow Long Baskets	17
<b>MediTray®   Insert Trays</b>	<b>18</b>
Endoscopy Trays	19
Insert Trays	21
<b>SteriTite®   Consumables</b>	<b>22</b>
Filters & Single Use Kits	23
Tamper Evident Seals	24
Indicators & Load Cards	25
<b>MediTray®   Case Trays</b>	<b>26</b>
Dual Use Case Trays	27
Insert Boxes	28
Wrappable Case Trays	29
Dental Cassettes	30
<b>MediTray®   Components &amp; Partitions</b>	<b>31</b>
Specialty Inserts	32
Backbone® Silicone Brackets	33
Backbone® Multi-Slot Silicone Brackets & Mat	34
Metal Brackets	35
Partitions & Posts	36
<b>Case Solutions®   Instrument Chemistries</b>	<b>37</b>
SuperNova® Cleaners	38
Cleaners, Solvents & Lubricant	39
Case Care Wipes & Cleaners	40
RO Water Systems	41
Pumps & Dispensers	42
Instrument Brushes & Brush Kits	43
Case Soil® Wash Indicator & Test Kits	44
<b>SteriTite®   Case Carts &amp; Shelving</b>	<b>45</b>
Closed Case Carts	46
Open Case Carts & Shelving	47
<b>CaseTrak360®   Software</b>	<b>48</b>
CaseTrak 360®	49
Additional Software Modules	50
<b>MediTray®   Customization</b>	<b>51</b>
Basic Sets	52
Web Sets®	53
Da Vinci® Instrument Sets	54
<b>Case Medical®   Custom Manufacturing</b>	<b>55</b>
Facilities Listing	56



## Company Overview

Case Medical is a U.S. manufacturer of medical reprocessing solutions that enable healthcare facilities to improve quality of care and patient outcomes, while managing and reducing costs. Our comprehensive approach ensures that all reprocessing systems work together seamlessly.

Simply put, our mission is to promote the health of your facility and your patients.

Case Medical, Inc., a vertically integrated manufacturer, is located in Bloomfield, NJ. We are ISO 13485 and ISO 27001 certified and FDA registered with FDA 510k clearance and CE mark for our products. The company has a partnership with the U.S. EPA under the Safer Choice program and participates in the Chemical Footprint Project to promote sustainability and safer chemistry.



Containers



Chemistry



Carts



Software





## What We Offer

Case Medical understands that our customers have many choices for medical device reprocessing solutions. Our customers tell us that they prefer to work with the Case Medical team and our solutions because of our integrated, systematic approach.

**360 Degree Approach:** Case Medical views everything that happens from the end of one surgery to the beginning of the next as one continuous process.

**Sustainability:** Case Medical is the only medical device reprocessing manufacturer ever recognized by the U.S. Environmental Protection Agency (U.S. EPA) as a Safer Choice Partner of the Year.

**Instrument Protection:** Case Medical containers and modular tray system ensure that surgical devices are safe, secure, and can be properly sterilized for patient care.

**Universal Compatibility:** Case Medical's container system has been cleared by the U.S. FDA for all current healthcare sterilization modalities and medical devices.

**Validated Products:** Case Medical products are validated through independent laboratories to perform as promised so you can be confident in our claims.

**Customization Program:** Case Medical will customize instrument sets to match surgeons' preferences in support of better patient care and more successful outcomes.

**UDI Integration:** 2D barcode integration enhances product functionality and links Case Medical products with our proprietary CaseTrak360 software.








## Our Vision

Case Medical provides only the highest quality, cost-effective products that are safe for patients, staff, and the environment. This enables healthcare facilities to achieve the holistic goals of the Cost, Quality and Outcomes (CQO) movement, established and championed by the Association for Health Care Resource & Materials Management (AHRMM).

- **Quality of Care** is confidence that every instrument is safe for patient use.
- **Cost Effective** medical reprocessing results in less waste, reduced inventory, and fewer OR delays.
- **Successful Outcomes** are achieved by using effective products that are safe for patients, staff, and the environment.
- **Successful Prevention** seeks to reduce the risks and costs associated with HAIs and hospital readmissions by choosing high quality, validated products.



# SteriTite®

## The Universal Sealed Container System

The SteriTite® container system is the ideal reusable packaging system. It is safer, more efficient and cost-effective than traditional wrap, and along with the MediTray® inserts, has been cleared by the U.S. FDA and has received the CE mark for use with all methods of sterilization commonly used in healthcare facilities.

All SteriTite containers are sold as a complete unit with lid, base and filter retention plates. SteriTite containers are competitively priced, virtually repair and replacement free, for a long useful life with potential for major cost savings.

- Full Size Containers
- Three-Quarter Size Containers
- Half Size Containers
- Long & Wide Containers
- Mini, Narrow & Narrow Long Containers
- FlashTite® Containers & Baskets

Note: For Biohazard transport, SteriTite sealed containers may be ordered without perforations in lid and base to isolate contaminated items.

## Full Size Containers



SteriTite® containers and components in various sizes

**Universal Compatibility:** SteriTite® containers are part of a universal system, meaning they are compatible with all medical devices and current sterilizers. Containers are available with either a perforated or solid base and feature color-coded ID plaques with 2D bar codes to make identification a snap. Use with corresponding baskets on page 13.

## Universal | Perforated Bottom with 4 filter retention plates



**Part ID**  
**Description** **SC04FG**  
4" High Full Size  
Universal Sealed  
Container  
**Dimensions**  
(LxWxH) 23.5"x11"x4"  
(60x28x10 cm)



**Part ID**  
**Description** **SC06FG**  
6" High Full Size  
Universal Sealed  
Container  
**Dimensions**  
(LxWxH) 23.5"x11"x6"  
(60x28x15 cm)



**Part ID**  
**Description** **SC05FG**  
5" High Full Size  
Universal Sealed  
Container  
**Dimensions**  
(LxWxH) 23.5"x11"x5 "  
(60x28x13 cm)



**Part ID**  
**Description** **SC08FG**  
8" High Full Size  
Universal Sealed  
Container  
**Dimensions**  
(LxWxH) 23.5"x11"x8"  
(60x28x20 cm)

## Solid Bottom | Solid Bottom with 2 filter retention plates



**Part ID**  
**Description** **SC04F**  
4" High Full Size  
Sealed Container  
**Dimensions**  
(LxWxH) 23.5"x11"x4"  
(60x28x10 cm)



**Part ID**  
**Description** **SC06F**  
6" High Full Size  
Sealed Container  
**Dimensions**  
(LxWxH) 23.5"x11"x6"  
(60x28x15 cm)



**Part ID**  
**Description** **SC05F**  
5" High Full Size  
Sealed Container  
**Dimensions**  
(LxWxH) 23.5"x11"x5 "  
(60x28x13 cm)



**Part ID**  
**Description** **SC08F**  
8" High Full Size  
Sealed Container  
**Dimensions**  
(LxWxH) 23.5"x11"x8"  
(60x28x20 cm)



## Three-Quarter Size Containers



SteriTite® containers are sized to accommodate loaner sets

**Drop-in Convenience:** DIN-sized SteriTite® containers eliminate the need to wrap custom graphic (loaner) trays. This saves time and reduces the risks, delays and costs associated with reprocessing punctured wrapped trays. Additional benefits include the ability to stack trays and process in steam as well as STERRAD® and other low temperature sterilizers. Use with corresponding baskets on page 14.

### Universal | Perforated Bottom with 3 filter retention plates



**Part ID** SC04QG  
**Description** 4" High 3/4 size Universal Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x4" (46x28x10 cm)



**Part ID** SC06QG  
**Description** 6" High 3/4 size Universal Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x6" (46x28x15 cm)



**Part ID** SC05QG  
**Description** 5" High 3/4 size Universal Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x5" (46x28x13 cm)



**Part ID** SC08QG  
**Description** 8" High 3/4 size Universal Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x8" (46x28x20 cm)

### Solid Bottom | Solid Bottom with 1 filter retention plate



**Part ID** SC04Q  
**Description** 4" High 3/4 size Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x4" (46x28x10 cm)



**Part ID** SC06Q  
**Description** 6" High 3/4 size Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x6" (46x28x15 cm)



**Part ID** SC05Q  
**Description** 5" High 3/4 size Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x5" (46x28x13 cm)



**Part ID** SC08Q  
**Description** 8" High 3/4 size Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x8" (46x28x20 cm)

## Half Size Containers



SteriTite® with barcode for tracking Instrument sets

**Secure Seal:** SteriTite® containers feature a knife edge fit between lid and base with unique clamping action to create a reliably secure seal. Containers also feature a barcode that enables the instrument set to be tracked to the patient level with our CaseTrak360® asset management software. Use with corresponding baskets on page 15.

## Universal | Perforated Bottom with 2 filter retention plates



**Part ID** SC04HG  
**Description** 4" High Half size Universal Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x4" (31x28x10 cm)



**Part ID** SC06HG  
**Description** 6" High Half size Universal Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x6" (31x28x15 cm)



**Part ID** SC05HG  
**Description** 5" High Half size Universal Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x5" (31x28x13 cm)



**Part ID** SC08HG  
**Description** 8" High Half size Universal Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x8" (31x28x20 cm)

## Solid Bottom | Solid Bottom with 1 filter retention plate



**Part ID** SC04H  
**Description** 4" High Half size Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x4" (31x28x10 cm)



**Part ID** SC06H  
**Description** 6" High Half size Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x6" (31x28x15 cm)



**Part ID** SC05H  
**Description** 5" High Half size Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x5" (31x28x13 cm)



**Part ID** SC08H  
**Description** 8" High Half size Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x8" (31x28x20 cm)



## Long & Wide Containers



D-ring handles enable aseptic lid removal for enhanced safety

**Safer Design:** SteriTite® containers and components are constructed of anodized aluminum, platinum-cured silicone, and passivated stainless steel hardware, making them less susceptible to performance issues that can result in improperly processed devices and OR delays. Use with corresponding baskets on page 16.

### Universal | Perforated Bottom with 4 filter retention plates



**Part ID** **SC04LG**  
**Description** 4" High Long Universal Sealed Container  
**Dimensions (LxWxH)** 27"x11"x4" (69x28x10 cm)



**Part ID** **SC08LG**  
**Description** 8" High Long Universal Sealed Container  
**Dimensions (LxWxH)** 27"x11"x8" (69x28x20 cm)



**Part ID** **SC06LG**  
**Description** 6" High Long Universal Sealed Container  
**Dimensions (LxWxH)** 27"x11"x6" (69x28x15 cm)



**Part ID** **SC05WG**  
**Description** 5" High Wide Universal Sealed Container  
**Dimensions (LxWxH)** 26.5"x16"x5" (67x41x13 cm)

### Solid Bottom | Solid Bottom with 2 filter retention plates



**Part ID** **SC04L**  
**Description** 4" High Long Sealed Container  
**Dimensions (LxWxH)** 27"x11"x4" (69x28x10 cm)



**Part ID** **SC08L**  
**Description** 8" High Long Sealed Container  
**Dimensions (LxWxH)** 27"x11"x8" (69x28x20 cm)



**Part ID** **SC06L**  
**Description** 6" High Long Sealed Container  
**Dimensions (LxWxH)** 27"x11"x6" (69x28x15 cm)



**Part ID** **SC05W**  
**Description** 5" High Wide Sealed Container  
**Dimensions (LxWxH)** 26.5"x16"x5" (67x41x13 cm)



## Mini, Narrow &amp; Narrow Long Containers



Visible case filter inside container lid

**Puncture Proof:** SteriTite® containers include smart features like visible filters and an offset vent pattern that is more resistant to strikethrough penetrations and less vulnerable to contamination and microbial intrusion. SteriTite Mini, Narrow and Narrow Long containers are ideal for scopes, robotic arms, bariatric and small instruments sets. Use with corresponding baskets on page 17.

### Universal | Perforated Bottom with 2 filter retention plates



**Part ID** SC02MG  
**Description** 2" High Mini Universal Sealed Container  
**Dimensions (LxWxH)** 12"x6"x2" (30x15x5 cm)



**Part ID** SC03MG  
**Description** 3" High Mini Universal Sealed Container  
**Dimensions (LxWxH)** 12"x6"x3" (30x15x8 cm)



**Part ID** SC04MG  
**Description** 4" High Mini Universal Sealed Container  
**Dimensions (LxWxH)** 12"x6"x4" (30x15x10 cm)



**Part ID** SC02NG  
**Description** 2" High Narrow Universal Sealed Container  
**Dimensions (LxWxH)** 18.5"x6"x2" (47x15x5 cm)



**Part ID** SC03NG  
**Description** 3" High Narrow Universal Sealed Container  
**Dimensions (LxWxH)** 18.5"x6"x3" (47x15x8 cm)

### Solid Bottom | 1 filter retention plate

**Part ID** SC02M  
**Description** 2" High Mini Sealed Container  
**Dimensions (LxWxH)** 12"x6"x2" (30x15x5 cm)

**Part ID** SC03M  
**Description** 3" High Mini Sealed Container  
**Dimensions (LxWxH)** 12"x6"x3" (30x15x8 cm)

**Part ID** SC04M  
**Description** 4" High Mini Sealed Container  
**Dimensions (LxWxH)** 12"x6"x4" (30x15x10 cm)

**Part ID** SC02N  
**Description** 2" High Narrow Sealed Container  
**Dimensions (LxWxH)** 18.5"x6"x2" (47x15x5 cm)

**Part ID** SC03N  
**Description** 3" High Narrow Sealed Container  
**Dimensions (LxWxH)** 18.5"x6"x3" (47x15x8 cm)

## Introducing the Narrow Long container for single, long and small groups of instruments

### Universal | Perforated Bottom with 4 filter retention plates



**Part ID** SC04NLG  
**Description** 4" High Narrow Sealed Container  
**Dimensions (LxWxH)** 27"x6"x4" (68.5x15x10cm)

**Part ID** SC05NLG  
**Description** 5" High Narrow Sealed Container  
**Dimensions (LxWxH)** 27"x6"x5" (68.5x15x12.7cm)



**Part ID** SC04NL  
**Description** 4" High Narrow Sealed Container  
**Dimensions (LxWxH)** 27"x6"x4" (68.5x15x10cm)

**Part ID** SC05NL  
**Description** 5" High Narrow Sealed Container  
**Dimensions (LxWxH)** 27"x6"x5" (68.5x15x12.7cm)

## FlashTite® Containers &amp; Baskets



FlashTite® valve plates for drier IUSS processing

**Rapid Sterilization:** FlashTite® filterless valve plates are a faster drying bacterial barrier alternative to our disposable filter systems making them ideal for Immediate Use Steam Sterilization (IUSS). FlashTite filterless valve plates are FDA cleared for both pre-vac and gravity steam sterilization. Use the corresponding basket to secure the contents.

## Containers

**Part ID** **SC06HGV**  
**Description** 6" High Half Size FlashTite Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x6" (31x28x15 cm)



**Part ID** **SC08HGV**  
**Description** 8" High Half Size FlashTite Sealed Container  
**Dimensions (LxWxH)** 12.3"x11"x8" (31x28x20 cm)



**Part ID** **SC06QGV**  
**Description** 6" High Mid Size FlashTite Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x6" (46x28x15 cm)



**Part ID** **SC08QGV**  
**Description** 8" High Mid Size FlashTite Sealed Container  
**Dimensions (LxWxH)** 18.1"x11"x8" (46x28x20 cm)



**Part ID** **SC06FGV**  
**Description** 6" High Full Size FlashTite Sealed Container  
**Dimensions (LxWxH)** 23.5"x11"x6" (60x28x15 cm)



**Part ID** **SC08FGV**  
**Description** 8" High Full Size FlashTite Sealed Container  
**Dimensions (LxWxH)** 23.5"x11"x8" (60x28x20 cm)



## Basket

BSKH15LH



BSKH35LH




BSKQ15LH



BSKQ35LH



BSKF15LH

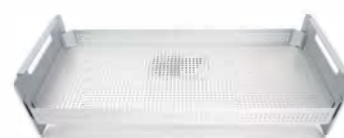

BSKF35LH







# MediTray®

## Instrument Baskets

The MediTray® family of instrument baskets fit within SteriTite® containers and feature a unique vent pattern that promotes fast drying. Additional system components including brackets, posts, partitions, and inserts, customize virtually any instrument set. Reconfigure as needed.

As with all SteriTite and MediTray products, aluminum components are anodized and stainless steel components are passivated for corrosion resistance.

- Full Size Baskets
- Three-Quarter Size Baskets
- Half Size Baskets
- Long & Wide Baskets
- Mini, Narrow & Narrow Long Baskets



## Full Size Baskets



MediTray® baskets in various sizes

**Versatile & Safe:** MediTray® instrument baskets are manufactured from lightweight anodized aluminum and feature design advantages including a highly perforated vent pattern for fast drying and flexible configuration. Instrument baskets are also offered in stainless steel to better withstand the damaging effects of caustic cleaners. Both are available in a variety of sizes. Use with corresponding containers on Page 6.

### Baskets | Full Size Aluminum



**Part ID**  
**Description**  
**Dimensions**  
(LxWxH)  
**BSKF01**  
1" High Full Basket  
21"x9.5"x0.9"  
(53.3x24.1x2.2 cm)



**Part ID**  
**Description**  
**Dimensions**  
(LxWxH)  
**BSKF02**  
2" High Full Basket  
21"x9.5"x1.8"  
(53.3x24.1x4.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(LxWxH)  
**BSKF03**  
3" High Full Basket  
21"x9.5"x2.8"  
(53.3x24.1x7 cm)



**Part ID**  
**Description**  
**Dimensions**  
(LxWxH)  
**BSKF04**  
4" High Full Basket  
21"x9.5"x3.7"  
(53.3x24.1x9.3 cm)



**Part ID**  
**Description**  
**Dimensions**  
(LxWxH)  
**BSKF05**  
5" High Full Basket  
21"x9.5"x4.7"  
(53.3x24.1x11.9 cm)

### Baskets | Full Size Stainless Steel

**BSKSF03**



SS Feet

**BSKSDF3**



SS End Panel



**Part ID**  
**Description**  
**Dimensions**  
(LxWxH)  
**BSKSF03/DF3**  
3" High Full  
SS Basket  
20.8"x9.8"x2.2"  
(52.8x24.8x5.5 cm)

"Stainless Steel  
baskets are  
available with  
attachable button  
feet or end panel  
identified with "D"

**BSKSF04**



SS Feet

**BSKSDF4**




SS End Panel



**Part ID**  
**Description**  
**Dimensions**  
(LxWxH)  
**BSKSF04/DF4**  
4" High Full  
SS Basket  
20.7"x9.4"x3.3"  
(52.5x23.8x8.3 cm)

## Three-Quarter Size Baskets

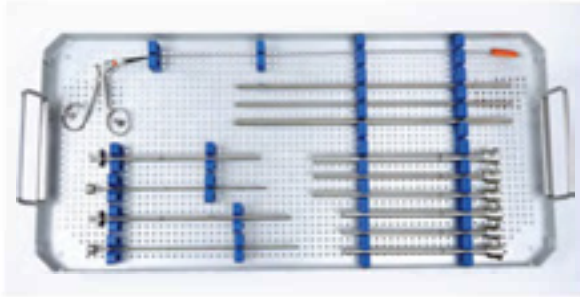

Customize your instrument sets by inserting components where needed

**Easily Configurable:** MediTray® Instrument baskets feature a unique vent pattern that not only promotes faster drying, it also provides a template for inserting Case Medical's proprietary system of brackets, posts, partitions, and inserts for flexible configuration. Available in anodized aluminum or passivated stainless steel for sterilization and decontamination of contents. Use with corresponding containers on page 7.

### Baskets | Three-Quarter Size Aluminum



**Part ID** BSKQ01  
**Description** 1" High 3/4 Basket  
**Dimensions (LxWxH)** 15.5"x9.5"x0.9" (39.3x24.1x2.2 cm)



**Part ID** BSKQ02  
**Description** 2" High 3/4 Basket  
**Dimensions (LxWxH)** 15.5"x9.5"x1.8" (39.3x24.1x4.5 cm)



**Part ID** BSKQ03  
**Description** 3" High 3/4 Basket  
**Dimensions (LxWxH)** 15.5"x9.5"x2.8" (39.3x24.1x7 cm)



**Part ID** BSKQ04  
**Description** 4" High 3/4 Basket  
**Dimensions (LxWxH)** 15.5"x9.5"x3.7" (39.3x24.1x9.3 cm)



**Part ID** BSKQ05  
**Description** 5" High 3/4 Basket  
**Dimensions (LxWxH)** 15.5"x9.5"x4.7" (39.3x24.1x11.9 cm)

### Baskets | Three-Quarter Size Stainless Steel

**BSKSQ03**



SS Feet

**BSKSDQ3**



SS End Panel

"Stainless Steel baskets are available with attachable button feet or end panel identified with "D"

**BSKSQ04**



SS Feet

**BSKSDQ4**



SS End Panel



**Part ID** BSKSQ03/DQ3  
**Description** 3" High 3/4 SS Basket with MediTray Feet  
**Dimensions (LxWxH)** 15.8"x9.8"x2.2" (40.1x24.8x5.5 cm)



**Part ID** BSKSQ04/DQ4  
**Description** 4" High 3/4 SS Basket with MediTray Feet  
**Dimensions (LxWxH)** 15.2"x9.4"x3.3" (38.6x23.8x8.3 cm)

## Half Size Baskets



Stackable baskets with convenient handles

**Aseptic Removal:** MediTray® Baskets secure contents during sterilization and handling. The stainless steel handles maintain their position and do not fall into loads allowing for safe, aseptic transfer of contents to the sterile field for patient procedures. While aluminum baskets dry quicker, stainless Steel baskets are also available for washing. Use with corresponding containers on page 8.

### Baskets | Half Size Aluminum



**Part ID** BSKH01  
**Description** 1" High Half Basket  
**Dimensions (LxWxH)** 9.5"x9.5"x0.9" (24.1x24.1x2.2 cm)



**Part ID** BSKH02  
**Description** 2" High Half Basket  
**Dimensions (LxWxH)** 9.5"x9.5"x1.8" (24.1x24.1x4.5 cm)



**Part ID** BSKH03  
**Description** 3" High Half Basket  
**Dimensions (LxWxH)** 9.5"x9.5"x2.8" (24.1x24.1x7 cm)



**Part ID** BSKH04  
**Description** 4" High Half Basket  
**Dimensions (LxWxH)** 9.5"x9.5"x3.9" (24.1x24.1x9.3 cm)



**Part ID** BSKH05  
**Description** 5" High Half Basket  
**Dimensions (LxWxH)** 9.5"x9.5"x4.5" (24.1x24.1x11.4 cm)

### Baskets | Half Size Stainless Steel

**BSKSH03**



SS Feet

**BSKSDH3**



SS End Panel

"Stainless Steel baskets are available with attachable button feet or end panel identified with "D"

**BSKSH04**



SS Feet

**BSKSDH4**



SS End Panel



**Part ID** BSKSH03/DH3  
**Description** 3" High Half SS Basket with MediTray Feet  
**Dimensions (LxWxH)** 9.3"x9.8"x2.2" (23.6x24.8x5.5 cm)



**Part ID** BSKSH04/DH4  
**Description** 4" High Half SS Basket with MediTray Feet  
**Dimensions (LxWxH)** 9.3"x9.4"x3.3" (23.6x23.8x8.3 cm)






## Long & Wide Baskets



Instrument baskets for every container size

**Large Device Organization and Protection:** MediTray® long and wide aluminum baskets are designed to accommodate longer and wider instruments without creating bends that can cause damage and prevent proper sterilization. Ideal for bariatric surgical devices, ring retractor systems, and flexible endoscopes. Available in stainless steel as a special order. Use with corresponding containers on page 9.

### Baskets | Long Size Aluminum

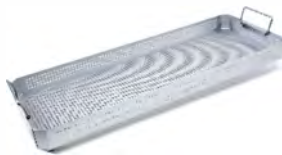

**Part ID** BSKL01  
**Description** 1" High Long Basket  
**Dimensions (LxWxH)** 25.3"x9.5"x0.9" (64.2x24.1x2.2 cm)



**Part ID** BSKL03  
**Description** 3" High Long Basket  
**Dimensions (LxWxH)** 25.3"x9.5"x2.8" (64.2x24.1x7 cm)



**Part ID** BSKL02  
**Description** 2" High Long Basket  
**Dimensions (LxWxH)** 25.3"x9.5"x1.8" (64.2x24.1x4.5 cm)



**Part ID** BSKL04  
**Description** 4" High Long Basket  
**Dimensions (LxWxH)** 25.3"x9.5"x3.7" (64.2x24.1x9.3 cm)

### Baskets | Wide Size Aluminum




**Part ID** BSKW02  
**Description** 2" High Wide Basket  
**Dimensions (LxWxH)** 24.7"x14.7"x1.8" (62.7x37.3x4.5 cm)



**Part ID** BSKW04  
**Description** 4" High Wide Basket  
**Dimensions (LxWxH)** 24.7"x14.7"x3.5" (62.7x37.3x8.8 cm)



## Mini, Narrow & Narrow Long Baskets

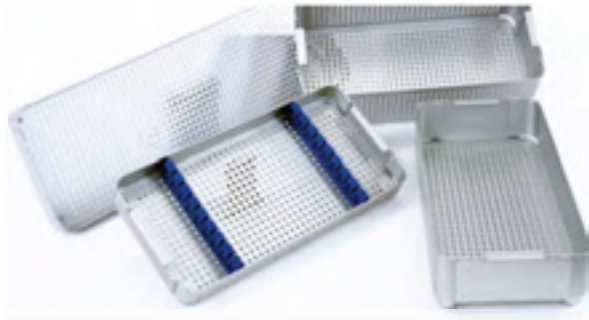

Unique vent pattern is a key MediTray® basket advantage

**Small Parts and Delicate Instruments:** MediTray® narrow baskets are sized to safely contain single devices, small instrument sets and delicate devices that could otherwise be lost or damaged in larger loads. The new Narrow Long baskets are ideal for scopes up to 25" long and for robotic arms. Use brackets, posts, partitions, and inserts to secure items in the baskets. Use with corresponding containers on page 10.

### Baskets | Mini Size Aluminum

	Part ID Description Dimensions (LxWxH)	<b>ST01</b> 1" High Mini Basket 11.5"x5.5"x0.9" (29.2x13.9x2.2 cm)
	Part ID Description Dimensions (LxWxH)	<b>ST02</b> 2" High Mini Basket 11.5"x5.5"x2" (29.2x13.9x5 cm)
	Part ID Description Dimensions (LxWxH)	<b>ST03</b> 3" High Mini Basket 11.5"x5.5"x3" (29.2x13.9x7.6 cm)

### Narrow Size Aluminum

	Part ID Description Dimensions (LxWxH)	<b>STN1</b> 1" High Narrow Basket 17.9"x5.5"x0.9" (45.4x13.9x2.2 cm)
	Part ID Description Dimensions (LxWxH)	<b>STN2</b> 2" High Narrow Basket 17.9"x5.5"x2" (45.4x13.9x5 cm)

Use Narrow Long baskets with SC04NL(G) and SC05NL(G) SteriTite Containers

### Baskets | Narrow Long.



Part ID	<b>STNL3</b>
Description	3" High Narrow Basket
Dimensions (LxWxH)	26.7"x5.7"x3" (67x14.4x7.6 cm)



Part ID	<b>STNL4</b>
Description	4" High Narrow Basket
Dimensions (LxWxH)	26.7"x5.7"x4" (67x14.4x10 cm)





# MediTray®

## Insert Trays

**S**ave time and space with stackable, interchangeable MediTray® insert trays. Reduce the need for duplicate containers and additional loads by providing a convenient way to secure instruments for a more efficient use of container space.

Trays fit within MediTray baskets which may be placed in compatible SteriTite® containers for sterilization, transport, storage, and aseptic removal of contents. Cleared for multi-modality sterilization and appropriate to wash in a bulk load.

— Endoscopy Tray Sets

— Insert Trays





## Endoscopy Tray Sets



Secure endoscopy instrument sets

**Secure instruments with “L” shaped handles:** The MediTray® Endo Tray System is designed for devices like graspers, kerrisons, and forceps. An interior bridge with scalloped, silicone brackets holds multiple instruments securely in place. This creates a more efficient use of space, maximizes exposure to sterilants, and promotes fast drying. The base tray is used to secure scopes, cameras, and other accessories. An insert tray (included) houses trocars, cannulas, and light cables.

### Insert Trays | Half Size Endo Trays



Part ID	<b>ENDO H1</b>
Description	Two-Tier System for use with 6" High Half Size Container
Dimensions (LxWxH)	9.8"x9.7"x5.7" (24.9x24.6x14.4 cm)

### Insert Trays | Three-Quarter Size Endo Trays



Part ID	<b>ENDO Q1</b>
Description	Four-Tier System for use with 8" High Three-Quarter Size Container
Dimensions (LxWxH)	16.1x9.7"x6.8" (40.9x24.6x17.2 cm)



Part ID	<b>ENDO Q1T234</b>
Description	Three-Tier System for use with 6" High Three-Quarter Size Container
Dimensions (LxWxH)	16.1x9.7"x5" (40.9x24.6x12.7 cm)



Part ID	<b>ENDO Q2</b>
Description	Two-Tier System for use with 4" High Three-Quarter Size Container
Dimensions (LxWxH)	16.1x9.7"x3.7" (40.9x24.6x9.4 cm)




Part ID	<b>ENDO Q1T456</b>
Description	Three-Tier System for use with 8" High Three-Quarter Size Container
Dimensions (LxWxH)	16.1"x9.7"x7" (40.9x24.6x17.8 cm)

Note: Use 3/4 size ENDO Q1 for an Arthroscopy set.

© Copyright Case Medical, Inc. Rev 2022



## Endoscopy Tray Sets

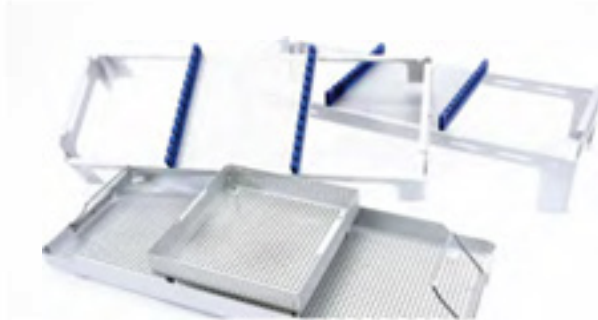

Containerize Lap Chole and Bariatric sets

**Multi-level stacking:** The MediTray® Endo Tray system reduces inventory and handling when multiple layers are used to keep everything needed for a procedure in one place at the same time. Place in a corresponding SteriTite® container.

### Insert Trays | Full Size Endo Trays



**Part ID** ENDOF1  
**Description** Four-Tier System for use with 8" High Full Size Container  
**Dimensions (LxWxH)** 21.6"x9.7"x7" (54.9x24.6x17.8 cm)



**Part ID** ENDOF1T456  
**Description** Three-Tier System for use with 8" High Full Size Container  
**Dimensions (LxWxH)** 21.6"x9.7"x7" (54.9x24.6x17.8 cm)



**Part ID** ENDOF2  
**Description** Two-Tier System for use with 4" High Full Size Container  
**Dimensions (LxWxH)** 21.6"x9.7"x3.7" (54.9x24.6x9.4 cm)




**Part ID** ENDOF1T234  
**Description** Three-Tier System for use with 6" High Full Size Container  
**Dimensions (LxWxH)** 21.6"x9.7"x5" (54.9x24.6x12.7 cm)

### Insert Trays | Long Size Endo Trays



**Part ID** ENDOL1  
**Description** Four-Tier System for use with 8" High Long Container  
**Dimensions (LxWxH)** 25.4x9.7"x6.7" (64.5x24.6x17 cm)

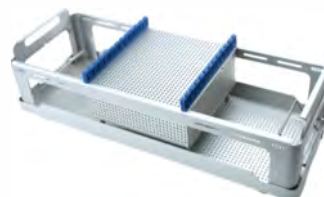

**Part ID** ENDOL1T456  
**Description** Three-Tier System for use with 8" High Long Container  
**Dimensions (LxWxH)** 25.4"x9.7"x6.7" (64.5x24.6x17 cm)

Note: Utilize the four tier ENDOF1 in a full size SteriTite container for Lap Chole when numerous instruments are required for a procedure.

© Copyright Case Medical, Inc. Rev 2022

## Insert Trays



Containerized tray conveniently lifts out as a single unit

**Less Risk of Damage:** MediTray® insert trays reduce the risk of damage to delicate and critical devices by providing dedicated containment that separates and secures these items during processing. Use with respective Meditray baskets. Posts are available to elevate and secure insert trays to baskets below.



Part ID **IMT01**  
Dimensions  
(LxWxH) 9.3"x6.2"x0.98"  
(23.6x15.7x2.4 cm)



Part ID **IHT01**  
Dimensions  
(LxWxH) 8.5"x9.3"x0.98"  
(21.5x23.6x2.4 cm)



Part ID **IQT01**  
Dimensions  
(LxWxH) 14.5"x9.3"x0.98"  
(36.8x23.6x2.4 cm)



Part ID **STX01**  
Dimensions  
(LxWxH) 17.2"x8.5"x0.98"  
(43.6x21.5x2.4 cm)



Part ID **ITX015**  
Dimensions  
(LxWxH) 17.2"x8.5"x1.6"  
(43.6x21.5x4 cm)



Part ID **IFT02**  
Dimensions  
(LxWxH) 19.8"x9.3"x2"  
(50.2x23.6x5 cm)



# SteriTite®

## Single Use Consumables

SteriTite® consumables ensure that the containerized set is tamper proof, that the sterility of contents is maintained during transport and handling, and that sets that have been sterilized are identified. Use paper filters for rapid processing in steam sterilization cycles. Polypropylene filters are universal and required for H2O2 sterilization.

Order by item or in kits suitable for the sterilization method and container model. A kit of filters, seals and load cards ensures that the necessary consumables are there as needed.

- Filters & Single Use Kits
- Tamper Evident Seals
- Indicators & Load Cards

## Filters &amp; Single Use Kits



*SteriTite® containers utilize single use filters*

**Protective Barriers:** SteriTite® sealed containers require single use filters to create a bacterial barrier. Use paper filters for rapid processing for terminal steam sterilization. Polypropylene filters are universal and required for H2O2 sterilization including use in STERRAD®, V-Pro®, Sterizone® and other hydrogen peroxide systems. Order a kit of filters, seals, and load cards and have everything you need for sterilization of SteriTite and other DIN-size containers. Square filters (9"x9") are also available from Case Medical.

## Consumables | Filters



**Part ID** SCF01  
**Description** Standard  
7.5" Diameter  
Cellulose Filters  
1000/pk



**Part ID** SCF02  
**Description** Standard  
7.5" Diameter  
Polypro Filters  
1000/pk



**Part ID** SCFM01  
**Description** Mini & Narrow  
10" x 4"  
Cellulose Filters  
1000/pk



**Part ID** SCFM02  
**Description** Mini & Narrow  
10" x 4"  
Polypro Filters  
1000/pk

## Consumables | Single Use Kits



**SCKIT1BP**  
Steam and EtO Consumable Kit  
Standard paper filters (SCF01)  
3 pk, 1000 ea/pk  
Blue seals (SCS01B) 1 pk, 1000 ea/pk  
Load cards (SCL01) 1 pk, 1000 ea/pk



**SCKIT1WN**  
H2O2 Consumable Kit  
Standard polypro filters (SCF02)  
3 pk, 1000 ea/pk  
White seals (SCS01W)  
1 pk, 1000 ea/pk  
Load cards (SCLH2023)  
1 pk, 1000 ea/pk



**SCKIT2BP**  
Steam and EtO Consumable Kit  
Mini paper filters (SCFM01)  
1 pk, 1000 ea/pk  
Blue seals (SCS01B) – 1 pk, 1000 ea/pk  
Load cards (SCL02) – 1 pk, 1000 ea/pk



**SCKIT2WN**  
H2O2 Consumable Kit  
Mini polypro filters (SCFM02)  
1 pk, 1000 ea/pk  
White seals (SCS01W)  
1 pk, 1000 ea/pk  
Load cards (SCLH2024)  
1 pk, 1000 ea/pk

**Note:** Use 7.5" diameter round filters for standard SteriTite and other DIN size containers. Use 10x4" rectangular filters for mini and narrow containers.

© Copyright Case Medical, Inc. Rev 2022



## Tamper Evident Seals



Color coded seals identify process

**Identification & Security:** SteriTite® container tamper evident seals provide assurance that sets are processed, identified, and tamper proof. Container seals are available for steam and low temperature sterilization cycles. Identifier seals are used to designate items requiring repair or sharpening. White seals may identify clean items. Red biohazard seals identify contaminated (used) items.

### Consumables | Tamper Evident Container Seals



Part ID  
Description

**SCS01B**  
Blue Container Seal  
for Steam & EtO  
Sterilization  
1000/pk



Part ID  
Description

**SCS01R**  
Red Container Seal\*  
for Steam & EtO  
Sterilization  
1000/pk



Part ID  
Description

**SCS01W**  
White Container  
Seal for H2O2  
Sterilization  
1000/pk

### Consumables | Identifier Seals



Part ID  
Description

**SCS01BH**  
Biohazard Seal  
1000/pk



Part ID  
Description

**SCS01W**  
White Seal  
for clean items  
1000/pk



Part ID  
Description

**SCS01RP**  
Repair Seal  
1000/pk



Part ID  
Description

**SCS01SH**  
Sharpen Seal  
1000/pk



## Indicators &amp; Load Cards



External load cards are sized to fit the respective container

**Visual Process Monitor:** SteriTite® Indicators and load cards are identification tags placed in dedicated slots on containers. The indicator changes colors when sets are sterilized. Steam and EtO indicators and load cards change from light to dark, while those for H2O2 change from red to orange (i.e. dark to light).


## Consumables | External Indicators



**Part ID Description** **SCL01**  
Load Cards  
(STEAM/EtO)  
Large for SteriTite  
Containers 1000/pk



**Part ID Description** **SCL02**  
Load Cards  
(STEAM/EtO)  
for Mini & Narrow  
SteriTite Containers  
1000/pk

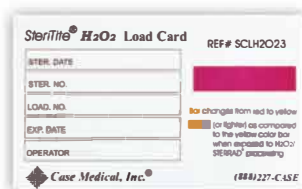

**Part ID Description** **SCLH2023**  
Load Cards  
(H2O2) for SteriTite  
Containers  
1000/PK



**Part ID Description** **SCLH2024**  
Load Cards  
(H2O2)  
for Mini & Narrow  
SteriTite Containers  
1000/pk

## Consumables | Internal Indicator



**Part ID Description** **SCI001**  
Dual Chemical  
Indicator for  
Steam and EtO  
1000/PK



# MediTray®

## Case Trays

The unique MediTray® system of case tray units with perforated lid and base may be used within MediTray baskets to provide secure containment of delicate items or specialty devices within a set. Use wrappable Case Trays when containers are not used.

Dual Use Case Trays provide the ultimate flexibility as they can be wrapped or placed within SteriTite® and other DIN sized sealed containers.

- Dual Use Case Trays
- Insert Boxes
- Wrappable Case Trays
- Dental Cassettes

## Dual Use Case Trays



*Wrap or drop-in convenience*

**Dual Purpose:** MediTray® Case Tray Inserts feature a lid, base, and locking mechanism. Use the provided screw-on feet when processed in SteriTite® and other DIN-sized sealed containers. Remove the feet to avoid tears when wrapped. BackBone® silicone brackets for securing rigid scopes are included with Scope Care Cases.



**Part ID** | **STCM2**  
**Description** | Mini Scope Care Case Tray w/ Silicone Brackets  
**Dimensions (LxWxH)** | 11"x5.5"x2" (28x14x5 cm)



**Part ID** | **STCN2**  
**Description** | Narrow Scope Care Case Tray w/ Silicone Brackets  
**Dimensions (LxWxH)** | 18"x5.5"x2" (46x14x5 cm)



**Part ID** | **HCI03**  
**Description** | Half Size Case Tray  
**Dimensions (LxWxH)** | 9.5"x9.5"x3" (24x24x7.6 cm)



**Part ID** | **QCI03**  
**Description** | Three-quarter Size Case Tray  
**Dimensions (LxWxH)** | 15.5"x9.5"x3" (39x24x7.6 cm)



**Part ID** | **FCI03**  
**Description** | Full Size Case Tray  
**Dimensions (LxWxH)** | 21"x9.5"x3" (53x24x7.6 cm)



**Part ID** | **LCI03**  
**Description** | Long Case Tray  
**Dimensions (LxWxH)** | 25"x9.5"x3" (63.5x24x7.6 cm)



**Part ID** | **WCI03**  
**Description** | Wide Case Tray  
**Dimensions (LxWxH)** | 24"x14"x4.5" (61x35x11.5 cm)





## Insert Boxes



Available in various sizes

**Grouping and Containment:** MediTray® insert boxes provide a secure, size-appropriate container that enables small components and accessories to be processed with larger instruments without risk of loss or damage. Insert boxes provide improved protection and sterilant penetration, and are a reusable alternative to peel packs which cannot be placed within a container or wrapped set.

	Part ID Dimensions (LxWxH)	<b>APC01</b> 5"x1.8"x1.0" (13x4.6x2.5 cm)		Part ID Dimensions (LxWxH)	<b>EC00</b> 10.3"x3.3"x0.8" (26.2x8.4x2 cm)
	Part ID Dimensions (LxWxH)	<b>APC02</b> 3.0"x3.0"x1.0" (7.6x7.6x2.5 cm)		Part ID Dimensions (LxWxH)	<b>EC01</b> 8.1"x4.5"x0.8" (20.6x11.4x2 cm)
	Part ID Dimensions (LxWxH)	<b>APC03</b> 5.1"x1.8"x2.0" (13x4.6x5.1 cm)		Part ID Dimensions (LxWxH)	<b>EC02</b> 9.8"x8.1"x0.8" (24.9x20.6x2 cm)
	Part ID Dimensions (LxWxH)	<b>APC04</b> 3.0"x3.0"x2.0" (7.6x7.6x5.1 cm)		Part ID Dimensions (LxWxH)	<b>EC04</b> 12.5"x9.8"x1.3" (31.8x24.9x3.3 cm)
	Part ID Dimensions (LxWxH)	<b>APC05</b> 5"x5.5"x2.7" (12.7x14x6.9 cm)		Part ID Dimensions (LxWxH)	<b>EC05</b> 12.7"x8.7"x0.97" (32.3x22.1x2.5 cm)
	Part ID Dimensions (LxWxH)	<b>APC09</b> 2.5"x1.8"x0.9" (6.4x4.6x2.3 cm)		Part ID Dimensions (LxWxH)	<b>EC07</b> 6"x3"x0.75" (15.2x7.6x1.9 cm)
				Part ID Dimensions (LxWxH)	<b>EC08</b> 2.4"x7.3"x1.7" (6.1x18.5x4.3 cm)



## Wrappable Case Trays



A full range of case trays based on your needs

**DIN size wrapped cases:** MediTray® case trays may be combined with medical grade wrap to create a sterile enclosure for individual devices, instrument sets or odd sized instruments that will not fit into a sealed container. Selected wrappable case trays are featured below. More available upon request.

### ● Wrapped Case Trays | Standard Size



**FC02**

Part ID  
Dimensions  
(LxWxH)

#### **FC02**

22.1"x10.6"x2.4"  
(56.1x26.9x6.1 cm)

#### **FC04**

22.1"x10.6"x4.3"  
(56.1x26.9x10.9 cm)



**QC04**

Part ID  
Dimensions  
(LxWxH)

#### **QC02**

16.6"x10.6"x2.4"  
(42.2x26.9x6.1cm)

#### **QC04**

16.6"x10.6"x4.3"  
(42.2x26.9x10.9cm)



**HC02**

Part ID  
Dimensions  
(LxWxH)

#### **HC02**

10.9"x10.6"x2.4"  
(27.7x26.9x6.1 cm)

#### **HC04**

10.9"x10.6"x4.3"  
(27.7x26.9x10.9 cm)



**MC02**

Part ID  
Dimensions  
(LxWxH)

#### **MC01**

11.6"x6.7"x1.4"  
(29.5x17x3.6cm)

#### **MC02**

11.6"x6.7"x2.2  
(29.5x17x5.6cm)

### ● Wrapped Case Trays | Narrow Size



#### **ENCM**

16.3"x4.1"x1.8"  
(41.4x10.4x4.6cm)



#### **XC03**

31.5"x7.2"x3.2"  
(80x18.3x8.1cm)

#### **XC03W11**

31.5"x11"x3.2"  
(80x28x8.1cm)

## Dental Cassettes



Dental cassettes per procedure

**Customizable Dental Trays:** MediTray® dental cassettes accommodate instruments for a variety of procedures, allowing the devices within the set to be processed together, saving time, and improving organization. Made from the highest quality anodized aluminum, the cassettes may be wrapped or placed inside a SteriTite® sealed container. Ideal for high volume clinics and dental schools. Order trays per procedure or by component.

## Dental Cassettes | Procedure Trays



**Part ID** DC00Y  
**Description** Crown, Bridge and Suture Removal  
**Dimensions (LxWxH)** 8"x4.7"x1.2" (20.3x11.9x3 cm)



**Part ID** DC01YD  
**Description** Basic Exam, Perio Personal, and Biopsy  
**Dimensions (LxWxH)** 8"x4.7"x1.8" (20.3x11.9x4.6 cm)



**Part ID** DC02YD  
**Description** Dental Hygiene  
**Dimensions (LxWxH)** 8"x5.9"x1.8" (20.3x15x4.6 cm)



**Part ID** DCOR  
**Description** Rubber Dam Custom Set  
**Dimensions (LxWxH)** 10.2"x7.9"x1.2" (26x20x3 cm)



**Part ID** DC04YD  
**Description** Oral, Surgical and Restorative  
**DC4GP** Grad Perio Custom Set  
**Dimensions (LxWxH)** 11.9"x9.6"x1.8" (30.2x24.4x4.6 cm)



**Part ID** DC03YD  
**Description** Cast Prep, Amalgam and Composite  
**Dimensions (LxWxH)** 10.7"x8"x1.8" (27x20.3x4.6 cm)

## Dental Cassettes | Components




**Part ID** DC01  
**Dimensions (LxWxH)** 8"x4.7"x1.8" (20.3x11.9x4.6 cm)  
**DT01** 7.6"x4.1"x.75" (19.3x10.4x1.9 cm)



**Part ID** DC03  
**Dimensions (LxWxH)** 10.7"x8"x1.8" (27x20.3x4.6 cm)  
**DT02** 8.37"x7.5"x.81" (21x19x2.0 cm)



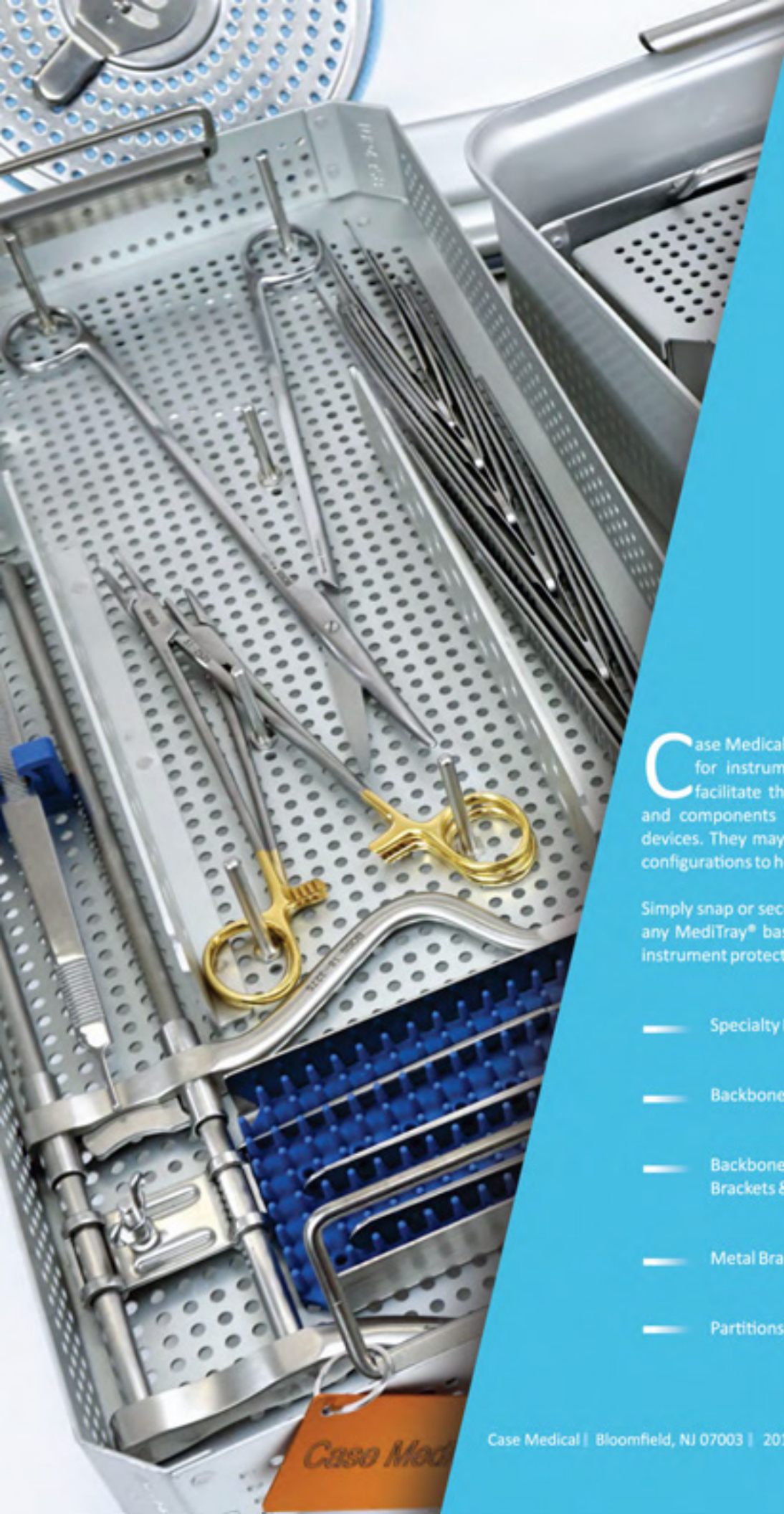

# MediTray®

## Components & Partitions

Case Medical manufactures numerous components for instrument protection, organization, and to facilitate the instrument count process. Brackets and components conform to the diameter of surgical devices. They may be used in endless combinations and configurations to hold instruments singly or in groups.

Simply snap or secure into the perforations at the base of any MediTray® basket, tray, or insert for the ultimate in instrument protection.

- Specialty Inserts
- Backbone® Silicone Brackets
- Backbone® Multi-Slot Silicone Brackets & Mats
- Metal Brackets
- Partitions & Posts

## Specialty Inserts



Assorted specialty inserts

**Versatile protection:** MediTray® accessories increase the value and versatility of the SteriTite® container system. Our wide range of options includes instrument stringers, screw caddies, and various racks for osteotomes, neuro dilators, and more. Consider our instrument stringers or instrument count rack for grouping and organizing ring-handled instruments. Use our peel pack rack to secure small packaged items and meet current standards.

### Specialty Inserts | Various options



**Part ID** APSRK01  
**Description** Screw Caddy  
**Dimensions (LxWxH)** 6.0"x4.8"x2.4"  
(15.2x12.2x6.1 cm)



**Part ID** APER01  
**Description** Rack for Ear Picks and Neuro Dilators  
**Dimensions (LxWxH)** 5.1"x5.6"x2.5"  
(13x14.2x6.4 cm)



**Part ID** AP0R1  
**Description** 10 Count Push Up Rack  
**Dimensions (LxWxH)** 9.8"x7.3"x0.9"  
(24.9x18.5x2.3 cm)

**AP0R2**  
20 Count Push Up Rack  
9.8"x7.3"x1.0"  
(24.9x18.5x2.5 cm)



**Part ID** APCR1  
**Description** Instrument Count Rack  
**Dimensions (Length)** 7.7"x9.8"x4.0"  
(19.5x24.9x10.2 cm)




**Part ID** INST06  
**Description** Instrument Stringer  
**Dimensions (Length)** 6" (15.2 cm)



**INST09**  
Instrument Stringer  
9" (22.9cm)



**APPW2**  
Peel Pack Rack  
13.5"x8.2"x1.0"  
(34.3x20.8x2.5cm)

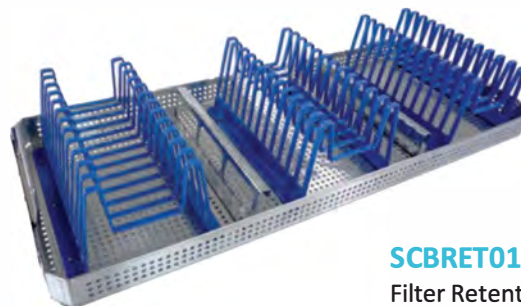

**SCBRET01**  
Filter Retention Plate  
24.9"x9.7"x4"  
(63.2x24.6x10cm)



## BackBone® Silicone Brackets



BackBone® brackets keep instrumentation in place

**Maximum Protection:** BackBone® silicone brackets organize and protect individual instruments during sterilization and transport. “Sure grip” features serve as shock absorbers with a patented inner spine and cushioning tabs to protect delicate devices that need to be held securely in place. Available in numerous sizes to secure devices from 1mm to 28mm diameter. Simply place and snap into the MediTray® perforation.

## Silicone Brackets | Single Slot



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BMS01**  
Silicone Small Bumper  
1.7"x0.4"  
(4.3x1.0 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BMS07**  
Single Slot Micro Brkt.  
5mm Diameter  
1.1"x1.0"  
(2.8x2.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BMS02**  
Single Slot Brkt.  
1-3mm Diameter  
1.1"x1.0"  
(2.8x2.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BES171**  
Single Slot Brkt.  
5-10mm Diameter  
1.1"x1.0"  
(2.8x2.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BES170**  
Single Slot Brkt.  
2-4mm Diameter  
1.1"x1.0"  
(2.8x2.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BMS12**  
Single Slot Brkt  
12mm Diameter  
1.5"x1.0"  
(3.8x2.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BES001**  
Single Slot Brkt.  
18mm Diameter  
1.8"x1.4"  
(4.5x3.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BES002**  
Single Slot Brkt.  
22mm Diameter  
1.9"x1.4"  
(4.8x3.5 cm)



**Part ID**  
**Description**  
**Dimensions**  
(WxH)

**BES003**  
Single Slot Brkt  
28mm Diameter  
2.2"x1.6" (5.5x4.0 cm)




## BackBone® Multi-Slot Silicone Brackets &amp; Mat



Non-slip, highly perforated MediTray® mat

**Secure Multiple Devices:** BackBone® silicone brackets and mats support and protect several instruments at a time within a MediTray® basket, tray or insert. BackBone scalloped brackets contain multiple slots to protect delicate devices in all sterilization modalities. Silicone mats provide an easy to use protective, non-slip layer with patented features that promote steam penetration during sterilization.

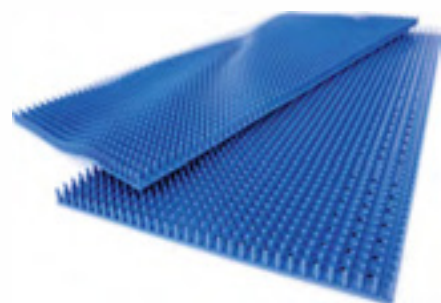

**Part ID** **MAT042**  
**Description** Full size silicone mats  
**Dimensions** 2 per pack 20.3"x9.5"  
 (LxWxH) (51.5x24.1cm)



**Part ID** **BMS10**  
**Description** 6 Slot Micro Brkt.  
**Dimensions** 5-7 mm Diameter  
 3.7"x.75"  
 (WxH) (9.4x1.9 cm)



**Part ID** **BMS03**  
**Description** Scalloped Brkt.  
**Dimensions** 5mm Diameter  
 (8 slots) 5.5"x1.0"  
 (WxH) (13.9x2.5 cm)



**Part ID** **BMS04**  
**Description** Scalloped Brkt.  
**Dimensions** 1-3mm Diameter  
 (8 slots) 5.5"x1.0"  
 (WxH) (13.9x2.5 cm)



**Part ID** **BMS05**  
**Description** Scalloped Brkt.  
**Dimensions** 5-10mm Diameter  
 (11 slots) 9.4"x1.0"  
 (WxH) (23.8x2.5 cm)



**Part ID** **BMS06**  
**Description** Scalloped Brkt.  
**Dimensions** 5-10mm Diameter  
 (12 slots) 9.4"x1.0"  
 (WxH) (23.8x2.5 cm)



**Part ID** **BMS08**  
**Description** Scalloped Brkt.  
**Dimensions** 2-4mm Diameter  
 (15 slots) 9.4"x1.0"  
 (WxH) (23.8x2.5 cm)

Note: BMS06 and BMS05 can be used separately in all MediTray® cases or as offset brackets in multi-level tray systems. Use BMS08 for multiple smaller diameter instruments.

© Copyright Case Medical, Inc. Rev 2022

## Metal Brackets



Metal brackets in a MediTray® basket

**Ideal for Heavy Instruments:** MediTray® metal brackets elevate and separate instruments to promote maximum sterilant exposure and facilitate the instrument count. Brackets are made of passivated stainless steel and designed to organize and protect medical devices. Available in numerous diameters and heights. Select from chart below.

## Single Brackets | Stainless Steel

Part ID	Description	Dimensions (WxH)
B01	STOP (Slot-less) Bracket	1.3"x1.8" (3.3x4.5 cm)
B02	6.3mm Slot Bracket	1.3"x1.8" (3.3x4.5 cm)
B03	9.5mm Slot Bracket	1.4"x1.8" (3.5x4.5 cm)
B04	12.7mm Slot Bracket	1.5"x1.8" (3.8x4.5 cm)
B05	15.8mm Slot Bracket	1.8"x1.8" (4.5x4.5 cm)
B06	19mm Slot Bracket	1.8"x1.8" (4.5x4.5 cm)
B07	25.4mm Slot Bracket	2.0"x1.8" (5.0x4.5 cm)
B08	31.7mm Slot Bracket	2.3"x1.8" (5.8x4.5 cm)

## Multi-Slot Brackets | Stainless Steel

Part ID	Description	Dimensions (WxH)
BM01	STOP (Slot-less)	2.9"x0.9" (7.3x2.2 cm)
BM02	8 Slot 5.5 mm Bracket use with Silicone Brackets BMS03 and BMS04	5.5"x0.9" (13.9x2.2 cm)
BM05	Partition	5.5"x0.9" (13.9x2.2 cm)
BM06 + BM07	5 Slot 8.5mm Brackets for grouping osteotomes and curettes	5.5"x0.9" (13.9x2.2 cm)
BM08	Partition	8.6"x0.9" (23.8x2.2 cm)
BM12	15 Slot 8 mm Bracket w/BMS08	9.4"x0.9" (23.8x2.2 cm)
BM15	11 Slot 14.7mm Bracket w/BMS05	9.4"x1.0" (23.8x2.5cm)
BM20	12 Slot 14.8mm Bracket w/BMS06	9"x1.1" (22.8x2.7cm)

## Partitions &amp; Posts



Partitions separate delicate from heavy devices

**Separate and Secure:** MediTray® partitions and posts separate and organize instruments within the tray. Partitions create compartments to separate instruments by type and to protect delicate instruments from damage. Posts provide separation between devices and may be used to stack ring-handled instruments in the open position. Attach to MediTray® base with threaded nuts available from Case Medical. Available in various sizes. See charts below.

## Partitions | Aluminum

Part ID	Dimensions (WxH)
P01	2.0"x1.8" (5.0x4.5 cm)
P02	3.0"x1.8" (7.6x4.5 cm)
P03	5.0"x1.8" (12.7x4.5 cm)
P04	9.4"x1.8" (23.8x4.5 cm)
P05	9.4"x3.7" (23.8x9.3 cm)
P06	9.4"x2.7" (23.8x6.6 cm)
P07	9.3"x1.5" (23.6x3.8 cm)
P08	2.0"x2.7" (5.0x7.1 cm)
P09	3.0"x2.7" (7.6x6.8 cm)
P10	5.0"x2.7" (12.7x6.8 cm)
P11	2.0"x3.6" (5.0x9.1 cm)
P12	3.0"x3.6" (7.6x9.1 cm)
P13	5.0"x3.6" (12.7x9.1 cm)



## Posts | Stainless Steel

Part ID	Dimensions (WxH)
POST01	0.9" (2.2 cm)
POST02	1.8" (4.5 cm)
POST03	2.7" (6.9 cm)
POST015	1.4" (3.5 cm)
POSTPLATE	3.8"x2" (9.6x5.0 cm)



Note: Aluminum Partitions are available in packages of 10.  
Stainless Steel Posts come in packages of 20.



# CASE Solutions® Alcohol Wipes

Pre-moistened with  
Ethanol & Isopropyl Alcohol

Single use alcohol wipe. Facilitates drying.  
Impregnated, non-woven material is  
pre-moistened with highly concentrated  
solution of isopropyl alcohol and ethanol in  
high purity H<sub>2</sub>O water.  
Quick, easy to use, non-streaking.  
Recommended for use in hospitals and  
healthcare facilities.

# CASE Solutions® Penta Wipes

Highly saturated with PentaZyme®

Penta Wipes are a pre-moistened, non-woven,  
disposable wipe. The single use wipe is  
impregnated with EPA preferred, biodegradable  
enzymatic cleaning solution.

Compatible with all medical devices.  
Easy to use, non-corrosive, pH neutral cleaner  
with a slight alcohol scent.

Recommended for use in hospital  
facilities

## Case Solutions® Instrument Chemistries

Case Medical instrument chemistries are validated by independent laboratories for outstanding cleaning effectiveness. Made from environmentally preferred ingredients, our cleaners and lubricant have been awarded the Safer Choice label by the U. S. Environmental Protection Agency.

With instrument chemistries for each step in the decontamination process, plus equipment and accessories to help you get the job done, call on Case Medical to improve your decontamination workflow.

- Supernova® Cleaners
- Case Solutions® Cleaners, Solvents & Lubricant
- RO Water Systems
- Pumps & Dispensers
- Instrument Brushes & Brush Kits
- Case Soil® Wash Indicator & Test Kits





## SuperNova® Cleaners



*Chemistries for every step*

**Super Concentrated for Super Efficiency:** SuperNova® concentrated liquid cleaners have been validated by independent laboratory testing to produce a 99.9% reduction in organic soil and 6 log reduction in bioburden, all with formulations that are safe for patients, staff, and the environment. SuperNova® cleaners are highly concentrated, so they cost less per use, and display the U.S. EPA Safer Choice label. Highly concentrated for use in automated washers. Follow with further sterilization or disinfection when processing surgical devices.



**Part ID** **CSNC01G**  
**Description** SuperNova® .1  
 Multi-Enzymatic Cleaner  
 Super Concentrate  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
**CSNC01W**  
 2x2.5 Gal (9.46 L) Btl per case



**Part ID** **CSNB01G**  
**Description** SuperNova® .1 NpHD  
 Neutral pH Detergent  
 Super Concentrate  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
**CSNB01W**  
 2x2.5 Gal (9.46 L) Btl per case



**Part ID** **CSNC25G**  
**Description** SuperNova® .25  
 Multi-Enzymatic  
 Cleaner Concentrate  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
**CSNC25W**  
 2x2.5 Gal (9.46 L) Btl per case



**Part ID** **CSNB25G**  
**Description** SuperNova® .25 NpHD  
 Neutral pH  
 Detergent Concentrate  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
**CSNB25W**  
 2x2.5 Gal (9.46 L) Btl per case



**Part ID** **CSNCD1G**  
**Description** SuperNova® Case Dry  
 Sheetting & Drying  
 Concentrate  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
**CSNCD1W**  
 2x2.5 Gal (9.46 L) Btl per case



**Part ID** **CSC011**  
**Description** PentaZyme®  
 Multi-Enzymatic  
 Cleaner Concentrate  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
 Ideal for Endoscopy Suite



**Part ID** **CSSTR011G**  
**Description** Schmutzoff® Descaler  
 Stainless Steel Cleaner  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
**CSSTR016W**  
 2x2.5 Gal (9.46 L) Btl per case



**Part ID** **CSNL01G**  
**Description** InstruCrema®  
 Lubricant Concentrate  
**Net Contents** 4x1 Gal (3.78 L) Btl per case  
**CSNL01W**  
 2x2.5 Gal (9.46 L) Btl per case





## Cleaners, Solvents & Lubricant



Pre-cleaners start the process

**Chemistries for Every Step:** Case Solutions® include an assortment of instrument chemistries, wipes, solvents, lubricant, and ready-to-use cleaners. Multi-enzymatic PentaPrep® or CasePrep®, our non-enzymatic pre-cleaner begin the process. Consider PentaZyme® with more active enzymes and our Bedside Kit for the endoscopy suite. Use PentaWipes® and Alcohol Wipes for surface cleaning, and Schmutzoff® for reconditioning stainless steel instruments, for removing residue from automated washers and autoclave walls. Apply CSR for ink and adhesive removing. Follow with subsequent cleaning and rinsing.



Part ID  
Description

### CSA011

PentaPrep®  
Multi-Enzymatic  
Pre-Soak Spray

Net Contents

12x8 oz (236 mL) Btl per case

### CSA015

12x24 oz (709 mL) Btl per case

### CSA012

4x1 Gal (3.78 L) Btl per case



Part ID  
Description

### CSP015

CasePrep™  
Non-Enzymatic  
Pre-Soak Spray

Net Contents

12x24 oz. (709 mL) Btl per case

### CSP011

12x8oz (236 mL) Btl per case

### CSP012

4x1 Gal (3.78 L) Btl per case



Part ID  
Description

### CSWE002

Penta Wipes®  
Multi-Enzymatic  
Surface Cleaner

Net Contents

Canister of 80 single use  
6 per case



Part ID  
Description

### CSWA002

Alcohol Wipes  
70% Denatured Alcohol  
Canister of 80 single use  
6 per case



Part ID  
Description

### CSSTR24

Schmutzoff® Descaler  
Stainless Steel Cleaner

Net Contents

12x24 oz. (709 mL)  
Btl per case



Part ID  
Description

### CSL014

InstruCreme®  
Lubricant

Net Contents

12x24 oz. (709 mL)  
Btl per case



### CSEKIT12

Endoscope Pre-Clean  
Bedside Kit 12x250mL bottles  
PentaZyme solution &  
12x3 pk. Penta Wipes per case



Part ID  
Description  
Net Contents

### CSR012

Ink & Adhesive Remover

12x8oz (236 mL) Btl per case

### CSR011

12x2oz. (59 mL) Btl per case





## Case Care®



Case Care® Wipes and Solutions for Cleaning and Sanitizing

**Safer Products for Personal Care and Infection Prevention:**

Formulated with only the safest ingredients applied to a single use wipe or a ready to use solution. From surface cleaning, cleansing, and sanitizing, Case Care® Wipes come in a handy, lightweight packet. Our infection prevention kit brings together the essentials for surface cleaning, sanitizing, and disinfection. Use the Case Care Cleanser or single use wipe when skin is soiled, before and after meals, and for personal care. Sustainable plant based cleanser with coconut, aloe, lavender and vitamin E.



**Part ID** CCW001  
**Description** Alcohol Wipes  
 Surface Disinfectant  
 70% Denatured Alcohol  
**Net Contents** 60 per pack (12 per case)  
**Dimensions** 8"x7" (20.3x17.7 cm)



**Part ID** CCW002  
**Description** Hand Sanitizer Wipes  
 70% Denatured Alcohol with  
 Aloe & Vitamin E  
**Net Contents** 60 per pack (12 per case)  
**Dimensions** 8"x7" (20.3x17.7 cm)



**Part ID** CCW003  
**Description** Surface Cleaner Wipes  
 Sustainable, pH Neutral  
**Net Contents** 60 per pack (12 per case)  
**Dimensions** 8"x7" (20.3x17.7 cm)



**Part ID** CCW004  
**Description** Skin Cleanser Wipes  
 Refresh hand and face  
**Net Contents** 60 per pack (12 per case)  
**Dimensions** 8"x7" (20.3x17.7 cm)



**Part ID** CSCOH24  
**Description** Surface Disinfectant  
 70% Denatured  
 Alcohol  
**Net Contents** 12x24 oz (709 mL)  
 Btl per case



**Part ID** CSHSAE03  
**Description** Hand Sanitizer  
 70% Denatured  
 Alcohol with  
 Aloe & Vitamin E  
**Net Contents** 12x24 oz (709 mL)  
 12x8oz (236 mL)  
 12x2oz (59 mL)  
 Btl per case



**Part ID** CSP015  
**Description** CasePrep™  
 Non-Enzymatic  
 Surface Cleaner  
**Net Contents** 12x24 oz (709 mL)  
 Btl per case



**Part ID** CSCCC24  
**Description** Case Care Cleanser  
 Hand & Body Wash  
**Net Contents** 12x24 oz (709 mL)  
 Btl per case

**CSKITIP02**

Infection Prevention Kit  
 12x 24 oz bottles with Applicator  
 Case Prep for Surface Cleaning,  
 70% Alcohol for Disinfecting  
 Case Care Hand Sanitizer for Infection Prevention



## RO Water Systems



RO system for optimal water purification

**Purer Water for Better Results:** Case Medical provides Reverse Osmosis (RO) Water Systems that efficiently remove contaminants including sediment, total dissolved solids (TDS), and microorganisms, that damage processing equipment and surgical devices. The result is critical water, required for the final rinse, which can improve the entire cleaning process and the steam quality in autoclaves, and enhances patient safety. Add an optional storage tank for consistent supply and a UV Light disinfection system.



**Part ID** CSROAT02  
**Description** AT Series RO system produces up to 500 gallons per day



**Part ID** CSROAT01  
**Description** AT Series RO system produces up to 1000 gallons per day



**Part ID** CSROBT01  
**Description** BT Series RO system produces up to 2,000 gallons per day



**Part ID** CSROCT01  
**Description** CT Series RO system produces up to 7,000 gallons per day



**Part ID** CSROTANK  
**Description** Storage Tank to collect the high-purity RO water for later use



**Part ID** CSROUV01  
**Description** UV Water Purification System destroys 99.9% of waterborne microorganisms



## Pumps & Dispensers



Accurate and economical dosing

**Accurate Dosing:** Case Medical offers automated pumps and dispensers to precisely control and record the amount of chemicals used in loads to produce the best possible results. They provide an electronic record of temperature and dosing for quality assurance and regulatory requirements. Automated dosing ensures that you always use the right concentration and get the biggest bang for your buck.



Pre-Soak Dosing System for manual cleaning

Part ID  
Description

**CSDMS01**  
Pre-Soak Dosing System with temperature sensor

**CSDMS02**  
Pre-Soak Dosing System with flow meter and temperature sensor

**CSDMS03**  
Pre-Soak Dosing System with USB port, flow meter and temperature sensor



Automated Foamer

Part ID  
Description

**CSFF**  
Fast Foam System for Pre-treatment



**SN01DSINK**  
Dispenses a chemical and water mixture into the sink  
Sliding a knob Foam Proportioner



Part ID  
Description

**CSCD03**  
3-product Liquid Chemical Dispenser with Remote Control and Washer Interface Module



## Instrument Brushes &amp; Brush Kits



Various sizes identified for use

**Cleaning Hard to Reach Areas:** Case Solutions® instrument cleaning brushes are perfect for lumened devices, crevices, and hard to reach surfaces. Nylon bristle brushes feature a proprietary blunted tip to minimize the risk of punctures and a twisted stainless wire handle for a rigid yet flexible grip. A color-coded bead identifies the diameter. Brushes are sold in packs of 12 or available in kits. Request the instrument brush rack to organize by type.

CSIBK1 GENERAL			CSIBK2 MICRO/EYE		
● Red	BR5F-06DZ	Brush 6" long, 5 french/1.67 mm	None	BR3F-06DZ	Brush 6" long, 3 french/1 mm
○ White	BR6F-08DZ	Brush 8" long, 6 french/2 mm	● Red	BR5F-06DZ	Brush 6" long, 5 french/1.67 mm
● Black	BR8F-08DZ	Brush 8" long, 8 french/2.6 mm	○ White	BR6F-08DZ	Brush 8" long, 6 french/2 mm
● Pink	BR10F-08DZ	Brush 8" long, 10 french, 3.3 mm	● Black	BR8F-08DZ	Brush 8" long, 8 french, 2.6 mm
● Orange	BR12F-08DZ	Brush 8" long, 12 french, 4 mm	● Pink	BR10F-08DZ	Brush 8" long, 10 french, 3.3 mm
○ White	BR03-12DZ	Twisted Wire, 12" long, 3 mm	● Orange	BR12F-08DZ	Brush 8" long, 12 french, 4 mm
● Green	BR05-16DZ	Twisted Wire, 16" long, 5 mm	● Red	BR02-12DZ	Twisted Wire, 12" long, 2 mm
● Yellow	BR07-16DZ	Twisted Wire, 16" long, 7 mm	○ White	BR03-12DZ	Twisted Wire, 12" long, 3 mm
● Pink	BR10-24DZ	Twisted Wire, 24" long, 10 mm			
CSIBK3 ENDOSCOPY					
● Red	BR02-12DZ	Twisted Wire 12" long, 2 mm	● Orange	BR2.5-24DZ	Twisted Wire 24" long, 2.5 mm
○ White	BR03-12DZ	Twisted Wire 12" long, 3 mm	● Green	BR05-24DZ	Twisted Wire 24" long, 5 mm
● Green	BR05-12DZ	Twisted Wire 12" long, 5 mm	● Pink	BR10-24DZ	Twisted Wire 24" long, 10 mm
● Green	BR05-16DZ	Twisted Wire 16" long, 5 mm	● Green	BR05-16DZ	Twisted Wire 16" long, 5 mm
● Pink	BR10-16DZ	Twisted Wire 16" long, 10 mm	● Pink	BR10-16DZ	Twisted Wire 16" long, 10 mm

Tooth brush style, endoscope, and double ended brushes are also available in packs of 12



Part ID | **CSIBK1**  
Description | General surgery



Part ID | **CSIBK2**  
Description | Micro eye



Part ID | **CSIBK3**  
Description | Endoscopy

Note: To clean your cannulated instruments, choose the correct size brush. The color of the brush identification bead is related to the diameter of the brush. The last two digits of the part number provide the length.



## Case Soil® Wash Indicator & Test Kits



Products to monitor the cleaning process

**Quality Assurance:** Case Soil® Wash Indicators and Case Solutions® Test Kits provide the ability to diagnose and address issues when results are not meeting expectations. Test Kits work at the front end of the process to assess water quality variables. Case Soil Indicators and Protein Detector provide information at the back end of the process to ensure that automated washers, ultrasonic cleaners, and detergents are working properly.

### Case Soil | Wash Indicator



**Part ID** **CSI001**  
**Description** Case Soil Wash Indicator  
100 per case



**CSI003**  
Case Soil Holder  
1 each



**CSI002**  
Protein Detector  
15 per case

### Solution Test Kits | Water Quality & Calibration



**Part ID** **CSTDs**  
**Description** TDS meter Indicates the total dissolved solids (TDS), pH and temperature of the water



**CSTEMPST**  
Single use temperature measure for automated washers



**Part ID** **CSCYL**  
**Description** Graduated cylinder kit to measure the volume of detergent drawn by automated washers



## SteriTite®

### Case Carts & Shelving

Case Medical manufactures open and closed case carts in single and double bay configurations, as well as triple bay storage racks. Designed for storage and transport, the universal size fits all sealed container systems and supplies. The stainless steel perforated shelves are fully adjustable and may be disassembled for cleaning.

The durable construction and materials withstand elevated temperatures and chemical agents. Storage racks are designed with the same durable, stainless steel materials and locking wheels.

- Closed Case Carts
- Open Case Carts & Shelving





## Closed Case Carts



*Stainless steel components for enduring quality*

**Secure Transportation and Storage:** The SteriTite® Case Cart system provides a safe and convenient way to transport and store surgical devices and supplies. Enhanced features include durable stainless steel construction, efficient use of space, adjustable perforated shelves, shock absorbing wheels with stainless steel casters, and a 2D bar code for tracking. Available in various sizes to meet the needs of the surgical specialty. Case carts with suffix D1 are extra deep and wide. Add a disposition monitor to identify clean from soiled (used) loads.



**CDS001** Identifies contents or use



**SCCC01**  
*Closed Single*



**SCCSD1**  
*Closed Deep Wide*



**SCCTD1**  
*Closed Deep High Wide*

Part ID	Description	Outer dimensions with Bumper DxWxH	Outer dimensions of Shell DxWxH	Number of Shelves	Depth and Width of Shelves
SCCC01	Closed Single	29"x22"x40"	26"x17"x33"	2	23"x12.6"
SCCC02	Closed Double	29"x38"x40"	26"x33"x33"	4	23"x12.6" (2 bays)
SCCSD1	Closed Deep Wide	34"x33"x40"	30"x27"x33"	2	28.3"x23.3"
SCCTS1	Closed Standard High	30"x33"x55"	26"x27"x48"	2	24"x23.3"
SCCTD1	Closed Deep High Wide	34"x33"x55"	30"x27"x48"	2	28.3"x23.3"



## Open Case Carts & Shelving



*Open Case Carts for storage and transport*

**Mobile Shelving and Transport:** SteriTite Open Case Carts & Shelving are designed for ease of use, and for storage and transport of containers and wrapped sets. Open case carts share the same enhanced features as closed case cart models, including adjustable perforated shelves to avoid tears in the wrappers, as well as shock absorbing wheels with stainless steel casters. No matter the size of the shelving unit it can be moved for transport and for environmental cleaning. Designed to house DIN-sized containers on individual shelves or to stack one on top of the other based on your requirements.



**SCCO02**

*Open Double Bay Cart*



**SCCO03H**

*Open Triple Bay Shelving*

Part ID	Description	Outer dimensions with Bumper DxWxH	Outer dimensions of Shell DxWxH	Number of Shelves	Depth and Width of Shelves
SCCO02	Open Double Bay Cart	27"x33"x44"	27"x33"x35"	4	26"x12.6"
SCCO03H	Open Triple Bay Shelving	27"x51"x65"	27"x51"x59"	9	26"x12.6"





## CaseTrak360® Software

CaseTrak360® software provides instrument tracking through every phase of reprocessing from decontamination to the operating room with full traceability between instrument sets and patients. It is all inclusive software at an affordable price, ideal for community based health care facilities as well as consolidated networks.

CaseTrak360 streamlines compliance with directives, improves communication, and shortens training periods and learning curves. Easily customizable to fit your work flow with touchscreen technology and fully cloud based. Utilize Loaner module, Tissue module, Endo module or Item module, as part of the CaseTrak360 software or acquire as individual modules. Use the UDI code on the container for instrument set tracking.

— CaseTrak360®

— Additional Software Modules

Note: For medical device companies requiring tracking and tracing of loaner sets, Case Medical offers its all inclusive software, Loaner Locate®.





## CaseTrak360®



*SteriTite® Container with 2D barcode and 2D set identifier*

**360° approach:** CaseTrak360® was developed to be the indispensable nerve center for high functioning instrument processing. Our extensive databases and proprietary algorithms consolidate all of the moving parts involved in sterile processing and the OR into informative dashboards, standardized reports and custom queries. CaseTrak360 uses an unique identification number to track instrument sets and specialty devices to the patient level for increased safety and infection prevention.




CaseTrak360 Main Menu provides navigation to each module and a hierarchy of use for each user.




CaseTrak360 Prep & Pack provides assembly instructions, direct link to oneSOURCE, photos, videos, count sheet, quality checks and alerts.



CaseTrak360 Dashboard View displays set information, i.e. location, turnaround, repair and missing items, and reprocessing instructions.

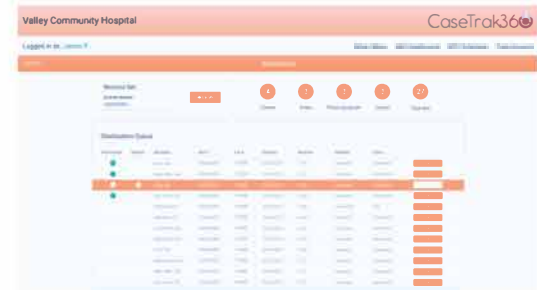

CaseTrak360 Sterilization module links equipment, monitoring, maintenance records, to a set for tracking and tracing.

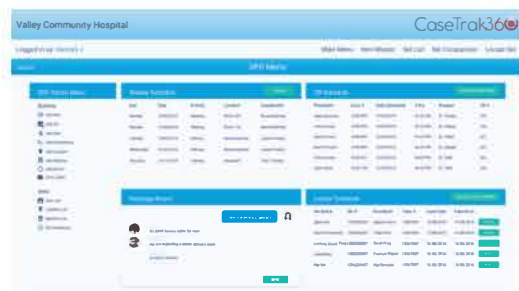

CaseTrak360 SPD Menu provides the schedule, actions and views to manage the department and administrate the functions.



CaseTrak360 OR module links a case to the OR schedule through staging, receiving, and return of items after a procedure

## Additional Software Modules



CaseTrak360® modules are easy to use, intuitive by design and comprehensive for the task.


**All that is required and more:** Loaner, Endo, Tissue and Item specific modules are additions to CaseTrak360® that can be used as stand alone solutions for the management and tracking of endoscopes, tissue and graft materials, and equipment assets. Designed to function individually they are a cost effective solution for specific tasks that don't require a software bundle. Now available, CaseTrak360 software for eye surgery centers to manage your process and your assets in real time.

## Software | Loaner Locate®



Loaner Locate®, finally, a software module designed to streamline and manage the loaner tray process from scheduling, receiving, use and return. Available to health care facilities as part of CaseTrak360 or as a stand alone software app.

## Software | Endo Module

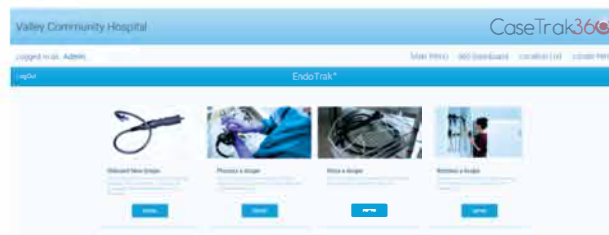

Endo Module tackles the workflow associated with endoscope re-processing, storage, and inventory management. As a stand alone module it is well suited to the needs of endoscopy clinics and larger facilities that handle endoscopes separate from the regular SPD.

## Software | Tissue Module



Tissue Module documents every step from receiving confirmation, storage and utilization of tissues and other biomaterials required for various surgical procedures including transplants and grafts.

## Software | Item Module



Item Module is the flexible software solution for tracking and managing individual items such as sterilization containers, equipment or medical devices. A bar coded item can be located simply and quickly in real time including any individual associated with its care and handling.



# MediTray®

## Customization

**M**editray® Customization sets the standard for the ultimate in instrument protection and organization. Our patented system of SteriTite® containers, trays and inserts provide flexible, modular layouts that are easy to use, reduce errors, and contribute to safer outcomes.

MediTray Basic Sets can be built for virtually any procedure. MediTray Web Sets® are pre-assembled, customized tray systems that provide the ultimate in instrument protection. Our representatives are trained to assist you in selecting a set for your needs.

- Basic Sets
- Web Sets®
- Da Vinci® Instrument Sets






SteriTite container (SC04QG) and basket (BSKQ03)

**Basic Customization:** Meditray® basic sets are a simple, cost-effective reusable option to containerize, sterilize, transport and secure surgical devices. Use the guide below to help you select the correct SteriTite® sealed container and basket(s) for a particular procedure tray.

Set Name	Container	Basket
AAA	SC05F(G)	BSKF04
Abdominal Hyst.	SC05F(G)	BSKF04
Amputation Set	SC04Q(G)	BSKQ03
Anterior Cervical	SC04Q(G)	BSKQ03
Appendix	SC06Q(G)	BSKQ02,BSKQ03
Arthroscope	SC03M(G)	ST02
Arthroscopy Set	SC06Q(G)	ENDOQ1T234
AV Fistula	SC04H(G)	BSKH03
AV Fistula Extras	SC04F(G)	BSKF03
Basic Heart Add	SC06F(G)	BSKF02, BSKF03
Basic Nasal	SC04Q(G)	BSKQ03
Basic Neuro	SC06F(G)	BSKF01, BSKF04
Basic Open Heart	SC08F(G)	BSKF03, BSKF04
Basic Plastic	SC04Q(G)	BSKQ03
Biopsy Set	SC04Q(G)	BSKQ03
Bone Mill	SC06H(G)	BSKH04
Bunionectomy	SC05H(G)	BSKH01, BSKH03
CABG	SC03M(G)	ST02
Carotid	SC04H(G)	BSKH03
Camera & Cord	SC04Q(G)	BSKQ03
Cataract Set	SC04Q(G)	BSKQ01, BSKQ02
Cervical Lam.	SC04Q(G)	BSKQ03
Chest Retractor	SC06Q(G)	BSKQ02,BSKQ03
Craniotomy	SC06Q(G)	BSKQ02,BSKQ03
C-Section Basic	SC04F(G)	BSKF03
Cysto	SC04Q(G)	ENDOQ2
D & C	SC04F(G)	BSKF03
Delivery	SC03M(G)	ST02
Dental Basic	SC04Q(G)	BSKQ01, BSKQ02
Ear Basic	SC06F(G)	BSKF02, BSKF03
Eye Tray	SC04Q(G)	BSKQ01,BSKQ02
Gastric Tray	SC06Q(G)	BSKQ01,BSKQ04
GU Lg. Instrument	SC04Q(G)	BSKQ03
GU Major	SC05Q(G)	BSKQ04

Set Name	Container	Basket
GU Minor	SC04H(G)	BSKH03
Hand & Foot	SC06Q(G)	BSKQ02,BSKQ03
Hand Set	SC06Q(G)	BSKQ02,BSKQ03
Head & Neck	SC06F(G)	BSKF02,BSKF03
Hemorrhoidectomy	SC05F(G)	BSKF04
Kidney Basic	SC06Q(G)	BSKQ02,BSKQ03
Laminectomy Retr.	SC06F(G)	BSKF02,BSKF03
Lami Neuro Basic	SC06F(G)	BSKF01,BSKF04
Laparoscopy tray	SC08F(G)	ENDOF1
Lap Chole	SC08F(G)	ENDOF1
Laparotomy	SC06F(G)	BSKF02,BSKF03
Large Plastic	SC04Q(G)	BSKQ03
Major Basic	SC06F(G)	BSKF02,BSKF03
Major Bone	SC04F(G)	BSKF01,BSKF02
Major Ortho	SC04F(G)	BSKF03
Major Trauma	SC05F(G)	BSKF04
Micro Plastic	SC03M(G)	ST02
Minor Basic Tray	SC04Q(G)	BSKQ03
Minor Ortho	SC04Q(G)	BSKQ03
Middle Ear	SC06F(G)	BSKF02, BSKF03
Myringotomy	SC04H(G)	BSKH03
Nasal Fracture	SC04Q(G)	BSKQ03
Neck & Head	SC04Q(G)	BSKQ03
Neuro Basic	SC06F(G)	BSKF02,BSKF03
Open Heart	SC08F(G)	BSKF03,BSKF04
Rectal/Hemorrhoid	SC04Q(G)	BSKQ03
Small Plastic	SC04H(G)	BSKH03
Sinus Endoscopy	SC06F(G)	BSKF02,BSKF03
T & A	SC04Q(G)	BSKQ03
Thoracotomy	SC05F(G)	BSKF04
Total Hip	SC05F(G)	BSKF04
Total Knee	SC05F(G)	BSKF04
Tracheostomy	SC04H(G)	BSKH03
Vaginal Hyst.	SC04F(G)	BSKF03

## Web Sets®

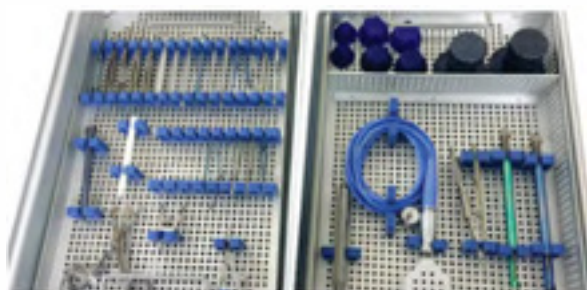


Cataract Set customized in SteriTite® sealed container

MediTray® Web Sets® are pre-assembled, customized tray systems that provide the ultimate in instrument protection and organization. Each tray system is designed to fit perfectly in its respective SteriTite® sealed container. Web Sets are configured for specific surgical specialties and specialty devices. With numerous options, our product specialists can assist you in selecting the perfect option.

Set Name	Container	Basket
Amputation Tray	SC04Q(G)	YRA1004D9C
Arthroscopy + Shaver	SC06F(G)	YRA1014A3A
Basic Orthopedic	SC04F(G)	YRA1027A4K
Biopsy Set	SC04Q(G)	YGA1001A
Camera Case	SC04H(G)	YRA1253A12H
Cataract Set	SC04Q(G)	YPA1002A3E
Cardio Tray	SC04F(G)	YNA1039A12E
Common Duct	SC04Q(G)	YGA1001A2W
Cysto Light & Cord Inst.	SC04FQ(G)	YSC1003A2W
D&C/L&D Set	SC04F(G)	YWA1002A2W
Ear Set	SC08F(G)	YQA1002A2W
Flexible Endoscope	SC04F(G)	YHA1024A9A
Hand and Foot	SC06F(G)	YRA1266A12N
Laminectomy	SC04Q(G)	YNA1013A4K
Large Ring Retractor Set	SC04F(G),SC05W(G)	BW01, BW04
Major Tray	SC08F(G)	YRA1234A12E
Minor Tray	SC05F(G)	YRA1233A9C
Myringotomy	SC04Q(G)	YQA1005B6N
Nasal Sinuscope	SC08F(G)	YQA1004A
Neuro Drill	SC06F(G)	YQA1082A9C
Neuro Set	SC06F(G)	YLA1000A3S
Ortho Tray	SC04F(G)	YRA1252A12E
Otoplasty	SC04Q(G)	YQA1101A13K
Sinuscopy Set	SC08Q(G)	YQA1007A2W
Sinus Lift Set	SC04H(G)	YTA1041A7S
Small Joint	SC04H(G)	YRA1008A2W
Strabismus	SC04Q(G)	YPA1057B13K
Tonsil and Adenoids	SC04F(G)	YQA1103A13L
Total Abd. Hysterectomy	SC08F(G)	YWA1001B13L
Trach Set	SC04F(G)	YQA1004A2W
Vascular Tray	SC06F(G)	YOA1004A5A



Cysto YSC1003A2W

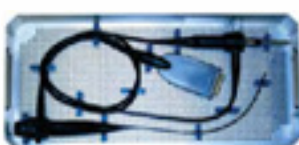Amputation Tray  
YRA1004D9CHand & Foot  
YRA1266A12NCamera Case  
YRA1253A12HFlexible Endoscope  
YHA1024A9AT & A Tray  
YQA1025C12LEar  
YQA1002A2WStrabismus  
YPA1057B13KOtoplasty  
YQA1101A13KMajor Basic  
YRA1234A12E






## Da Vinci® Instrument Sets




SteriTite container & custom insert for robotic arms.


**Developed with your input and validated with Intuitive Surgical:** Introducing a containerized option to sterilize, transport, and secure your valuable da Vinci robotic devices, instrumentation, and endoscope. No need to wrap these devices and be subject to tears, punctures, and subsequent reprocessing. Our SteriTite containers and custom trays provide the utmost device protection and assures that the device is ready when needed in the OR for the next procedure.

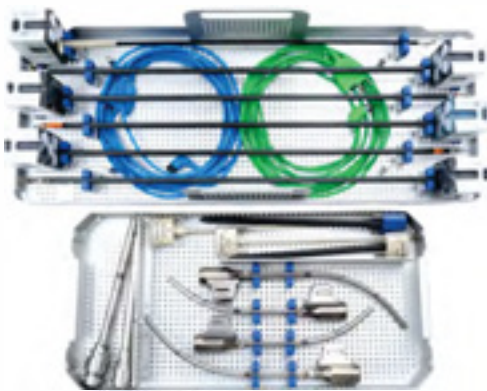

Da Vinci Xi Single Site  
YNB1068B16N + SC06L(G)

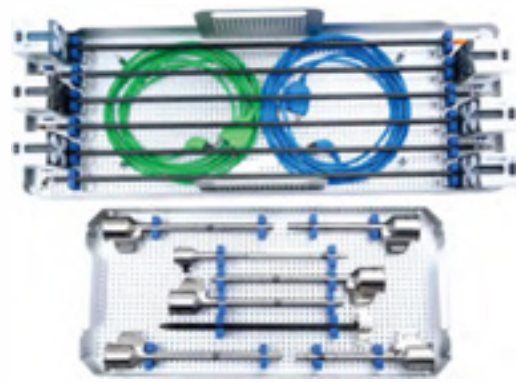

Da Vinci Xi Instrument Set  
YNB1067B16N + SC06L(G)

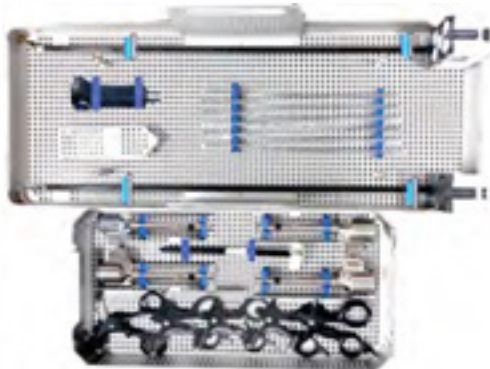

Custom Da Vinci Xi Instruments  
YNB1067B16N + SC06L(G)



Da Vinci Xi Endoscope Tray  
YNB1009O21H + SC05W(G)



Da Vinci Xi One arm tray  
YNB1177D21L + SC04NL(G)



Da Vinci Xi Three arms tray  
YNB1177E21N+ SC05NL(G)





# Case Medical®

Custom Manufacturing

Case Medical is an ISO certified, FDA registered medical device manufacturer. Our team of engineers utilize state-of art CAD/CAM technology, including 3-D modeling. All designs and products are produced in-house where all aspects are controlled and custom-tailored for exact compliance to customer's specifications and regulatory requirements.

Bringing all processes in-house, we've saved considerable time in providing you with the best quality at an affordable cost.

— Facilities Listing



## Facilities Listing



*Advanced technology for quality products*

**Vertically integrated manufacturer:** Case Medical® has invested heavily in the latest technology. Technologically advanced equipment include: hydraulic punch presses, CNC lasers, 450 Ton drawing press, complete machine shop, welding, and anodizing/passivation line. We are experienced in all phases of sheet metal manufacture, liquid injection molding, chemical formulation, and software development. The following facilities listing is to inform you of our capabilities.

### ENGINEERING:



State of the art 3-D CAD/CAM system with capability to accept all major CAD files.

Advanced Manufacturing Resource Planning (MRP) system.

Direct link from engineering to shop machinery.

### MANUFACTURING:

Amada Vipros King CNC turret press  
 Amada Vipros ACE CNC turret press  
 Amada Pulsar 1212 CNC lasers  
 CNC nylon coating machine  
 Hydraulico 450 Ton drawing press  
 Omera CNC trimming machine  
 Star SH12 CNC Swiss screw machine with automatic feed  
 Hardinge CNC lathe  
 Milltronics 3 axis machining center  
 Amada CNC press brakes with 3 axis back gauge  
 Wysong 5' shear  
 Mechanical press brakes  
 Haegar riveting/press machines  
 Empire bead blasting machines  
 Timesavers 37" wet surface sander  
 Falls edge deburring machine  
 Spot welders including aluminum spotwelding  
 NCG Sureweld TIG welder  
 Vibratory Tumbler/polishers  
 Ryobi cut off saw  
 Do All band saw  
 Floor lathe  
 Fray milling machine  
 Niagara circle cutter  
 Silk screening stations, image processor and drying oven  
 Laser Marking  
 Epson Laser Etching Machine  
 Engraving machine and chemical etch station  
 Automatic filling, labeling & capping line EE-440B  
 QC lab with spectrophotometer  
 Silicone Injection Molding Machine  
 3D Printer for custom prototypes  
 Anodizing and Passivating Line  
 Automated Wet Wipe Machine  
 CNC Controlled Filter Cutting Equipment

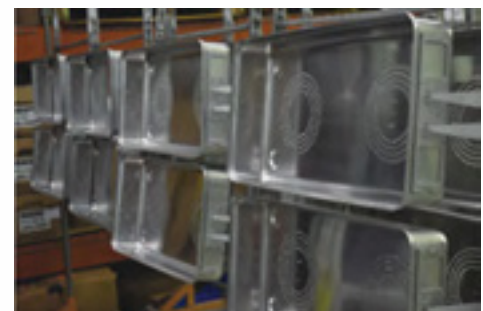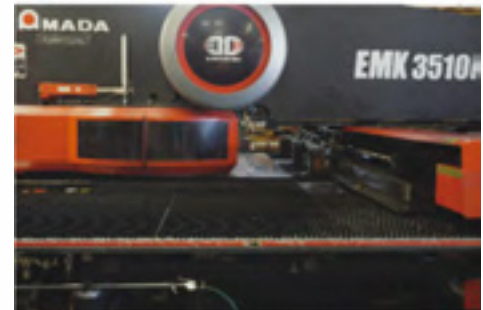



***Our mission is to promote the health  
of your facility and your patients***

Backbone, Case Soil, Case Solutions, CaseTrak360, FlashTite, InstruCrete, LoanerLocate, MediTray, PentaPrep, PentaWipes, PentaZyme, Schmutzoff, SteriTite, SuperNova are registered trademarks of Case Medical.

Da Vinci is a registered trademark of Intuitive Surgical.

SteriZone is a registered trademark of TSO3 Inc.

STERRAD is a registered trademark of Advanced Sterilization Products, Division of Ethicon US, LLC

V-PRO is a registered trademark of STERIS plc.

oneSOURCE is a trademark of Best Practice Professionals.

Case Medical, Inc.  
50 West Street,  
Bloomfield, NJ 07003  
U.S. and Canada 1-888-227-CASE  
Worldwide +1 201-313-1999  
[www.casemed.com](http://www.casemed.com)

© Case Medical, Inc. All rights reserved.



For more information contact  
[info@casemed.com](mailto:info@casemed.com)  
Phone: 888-227-CASE (2273)  
[www.casemed.com](http://www.casemed.com)



Scan QR code





For more information contact



**Case Medical**

50 West Street, Bloomfield, NJ 07003

info@casemed.com

Phone: 888-227-CASE (2273)

www.casemed.com



Scan QR code